



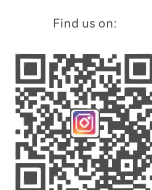

AP-A Dental Air Polisher  
More efficient scaling

Guilin Woodpecker Medical Instrument Co., Ltd.  
Information Industrial Park, Guilin National High-Tech  
Zone, Guilin, Guangxi, 541004 P.R. China  
Sales Dept.: +86-773-5873176  
E-mail: woodpecker@glwoodpecker.com  
Website: http://www.glwoodpecker.com  
ZHEUAC-012

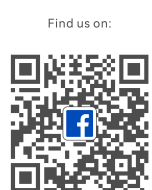

## AP-A

Sleek design

Cloud white  
appearance

Peerless elegance  
and practicality

Stand out  
where others  
fade into the masses


A unique Dental Air  
Polisher is in!



|                                                          |                                                    |
|----------------------------------------------------------|----------------------------------------------------|
| Main unit<br>1 unit                                      | Foot pedal<br>1 pc                                 |
| Power adapter (power cord included)<br>1 pc              | Subgingival air polishing handpiece (AP-2)<br>1 pc |
| Supragingival air polishing handpiece (AP-1)<br>1 pc     | Subgingival powder tank                            |
| Supragingival powder tank<br>1 pc                        | 25µm subgingival glycine powder<br>1 bottle        |
| 65µm supragingival sodium bicarbonate powder<br>1 bottle | Nozzle wrench<br>1 pc                              |
| Subgingival nozzle<br>10 pcs                             | 1400mL water bottle<br>1 pc                        |
| Sterilization box<br>2 pcs                               | Needles (Long & Short)<br>1 kit                    |
| 700mL conical water bottle<br>1 pc                       | Male and female connectors<br>1 pair               |
| Union tee<br>1 pc                                        | Air pipe connector<br>1 pc                         |



## Faster and more convenient scaling

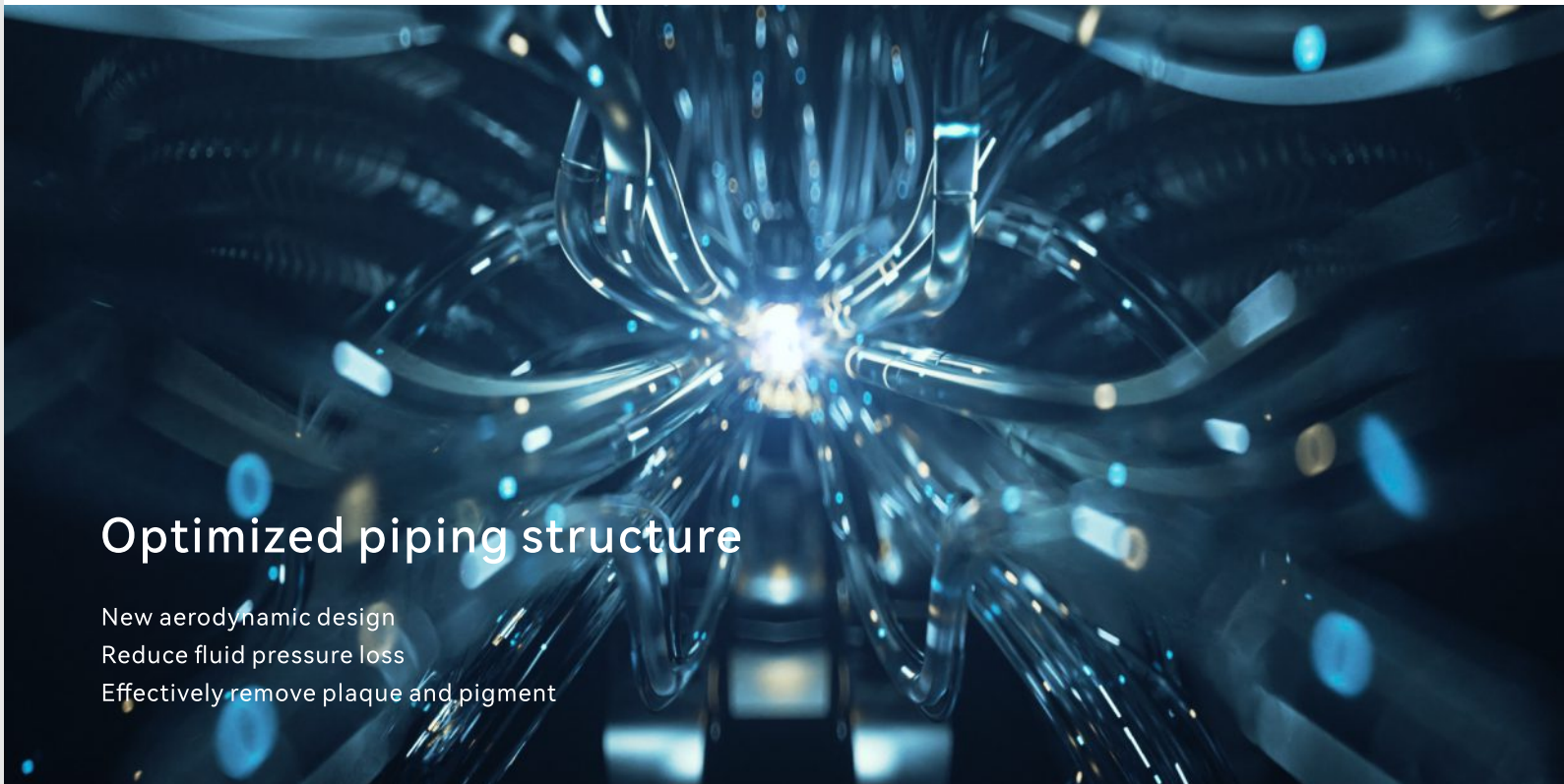

### Optimized piping structure

New aerodynamic design  
Reduce fluid pressure loss  
Effectively remove plaque and pigment



Before scaling



Plaque disclosing



After scaling

## Classic sodium bicarbonate powder

Average particle size:  
about 65 microns; hardness is less than that of  
enamel, doing no harm to teeth  
Contain natural plant-derived essence for a more  
pleasant taste  
The powder is waterproof to prevent the device from  
being blocked by wet powder

| Teeth    | Moh's hardness |
|----------|----------------|
| Enamel   | 4-5            |
| Dentin   | 2-2.5          |
| Cementum | 2-4            |

| Powder             | Moh's hardness |
|--------------------|----------------|
| Sodium bicarbonate | 2.5            |
| Glycine            | 2              |
| Calcium carbonate  | 3              |
| Aluminum hydroxide | 3              |

The promotion and sale of SUPER and SUPER PRO are not aimed at  
the following countries: France, Spain, Germany, Switzerland,  
Norway, Czech Republic, United Kingdom, Italy, Netherlands, Austria,  
Belgium, Denmark, Turkey, Sweden, Liechtenstein, Japan, and the  
United States.




## Two-in-one: Supragingival & subgingival mode

### Three-section handpiece design



Small-caliber nozzle  
More concentrated powder output  
More accurate scaling

Detachable handpiece,  
easy for negotiation if clog occurs  
Lightweight design for  
lower hand fatigue

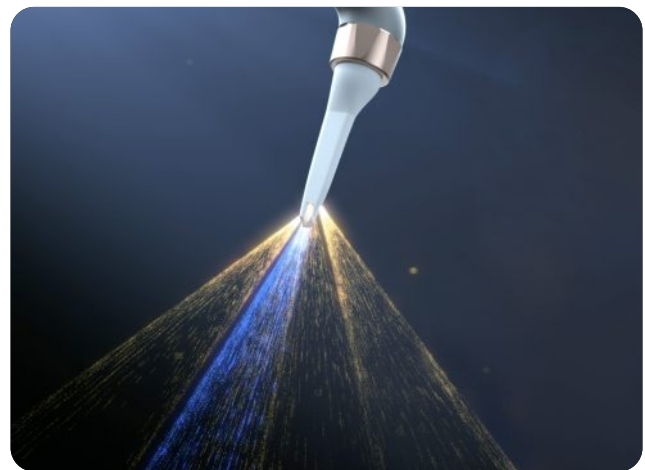

Subgingival nozzle made of medical polymers  
360° rotatable  
Bendable

Three-direction powder output  
One-direction water output  
Thorough plaque removal  
Great accessibility

## Ergonomic design

### Fully transparent funnel-shaped powder tank

Stable powder output  
360° rotatable, allowing real-time observation of  
powder flow in the tank

### Compact and lightweight

The main unit weighs only 2kg  
Space-saving base design

